



MotionPro DAS
User Manual
(For Windows™)

Software Release

1.05

Document Revision

September 2015

Products Information

<http://www.idtvision.com>

North America

1202 E Park Ave
TALLAHASSEE FL 32301
United States of America
P: (+1) (850) 222-5939
F: (+1) (850) 222-4591
llourenco@idtvision.com

Europe

via Pennella, 94
I-38057 - Pergine Valsugana (TN)
Italy
P: (+39) 0461- 532112
F: (+39) 0461- 532104
pgallorosso@idtvision.com
Eekhoornstraat, 22
B-3920 - Lommel
Belgium
P: (+32) 11- 551065
F: (+32) 11- 554766
amarinelli@idtvision.com

Copyright © Integrated Design Tools, Inc.

The information in this manual is for information purposes only and is subject to change without notice. Integrated Design Tools, Inc. makes no warranty of any kind with regards to the information contained in this manual, including but not limited to implied warranties of merchantability and fitness for a particular purpose. Integrated Design Tools, Inc. shall not be liable for errors contained herein nor for incidental or consequential damages from the furnishing of this information. No part of this manual may be copied, reproduced, recorded, transmitted or translated without the express written permission of Integrated Design Tools, Inc.

Table of Contents

1. TERMS AND CONDITIONS.....	5
2. PRECAUTIONS.....	5
2.1. Cables.....	5
3. SYSTEM OVERVIEW.....	6
3.1. Introduction to the MotionPro Data Acquisition System.....	6
3.2. System components.....	8
3.3. Note on operating systems.....	8
3.4. Software Development Kit.....	9
4. INSTALLING THE MOTIONPRO DAS.....	10
4.1. Minimum computer requirements.....	10
4.2. Package contents.....	10
4.3. Software Installation.....	11
4.4. Hardware Installation.....	11
4.5. Data Acquisition Back panel.....	12
5. DATA ACQUISITION STAND-ALONE PROGRAM.....	13
5.1. Open data acquisition device.....	13
5.2. Data Acquisition Menu structure.....	15
5.3. FILE Menu.....	16
5.3.1. Open a data file.....	17
5.3.2. Save acquired data.....	18
5.4. PROCESS Menu.....	19
5.4.1. Filters and Math operations.....	20
5.4.2. FFT and Power Spectrum.....	21
5.5. DEVICE Menu.....	22
5.5.1. Run Control.....	22
5.5.2. Vertical controls.....	23
5.5.3. Horizontal controls.....	24
5.5.4. Input trigger.....	25
5.5.5. Input channels.....	26
5.5.6. Configuring input channels.....	27
5.5.7. Output channels.....	28
5.6. View Menu.....	30
5.6.1. Average and Standard Deviation.....	30
5.6.2. Cursors.....	31
5.7. TOOLS menu.....	32
5.7.1. Waveform Editor.....	32
5.7.2. Options.....	32
5.7.3. Language selection.....	33
5.8. WINDOW Menu.....	33
5.9. HELP Menu.....	33
6. WAVEFORM EDITOR.....	34
6.1. Overview.....	34
6.2. The waveform editor menu structure.....	37
6.3. File Menu.....	38
6.4. Create a new waveform.....	38
6.5. Open a waveform file.....	39
6.6. Save a waveform.....	39
6.7. Edit Menu.....	40

6.8. Add components to a waveform.....	41
6.9. Insert primitives into a component.....	42
6.10. Edit a component.....	43
6.11. Delete a component.....	44
6.12. View Menu.....	45
6.13. Tools Menu.....	45
6.13.1. Options.....	45
6.13.2. Language.....	45
6.14. Help Menu.....	46
7. APPENDIX A – PRODUCT SPECIFICATIONS.....	47
7.1. Specifications.....	47

1. Terms and Conditions

For more detailed information, see the “**Terms and Conditions**” as stated in the IDT web site (www.idtvision.com/solutions/manuals.php).

2. Precautions

2.1. Cables

Ensure that all cable connections are properly secured and that there is not excessive strain on the cabling.

3. System Overview

3.1. Introduction to the MotionPro Data Acquisition System



The MotionPro Data Acquisition System (DAS) is a multifunction data acquisition module for the USB 2.0 bus. The module can perform simultaneous operation of analog input and analog output. The key hardware functions are listed below:

- 16 analog input channels with programmable gains.
- 4 analog output channels for waveform generation.
- Internal or external clock source.
- Trigger operations using a software command, an analog threshold value, or an external digital input value as the trigger event.
- 500 V galvanic isolation barriers that prevents ground loops to maximize analog signal integrity and protect your computer.

Analog input subsystem

- 16-bit A/D converter.
- Throughput rate up to 500 kHz.
- 16 single ended analog input channels.
- Programmable gain of 1, 2, 4 or 8 provides input ranges of ± 10 , ± 5 , ± 2.5 and ± 1.25 V.
- 1024-location channel-gain list. You can cycle through the channel-gain list by using continuous scan mode or triggered scan mode. The maximum sampling rate when using the channel-gain list is 500 kHz.

Analog output subsystem

- 16-bit D/A converter.
- Output rate up to 500 k samples/s.
- Output range of ± 10 V.

- The DACs are deglitched to prevent noise from interfering with the output signal.
- 1024 location output channel list. You can cycle through the output channel list using continuous output mode or waveform mode. For waveform generation mode, you can simultaneously update all four DACs at 500 kHz per channel; for continuous output mode, you can simultaneously update all four DACs at 250 kHz per channel.
- External or internal clock source.
- Trigger operations using a software command, an analog threshold value, or an external digital input value as the trigger event.
- A 500 V galvanic isolation barrier that prevents ground loops to maximize analog signal integrity and protect your computer.

3.2. System components

The MotionPro Data Acquisition System components are listed below.

- **Data Acquisition System module:** The MotionPro DAS module provides 16 analog inputs, 4 analog outputs, and 16 digital I/O with internal or external clock source and triggers.
- **Digital Interface:** The USB 2.0 (480 Mb/sec) cable provides data and control signals to and from the module. It connects the module to any USB 2.0 port on your personal computer or laptop.
- **Power Source:** The power supply provides external power to the device (5 V).
- **Data Acquisition Software and SDK:** it operates in Windows 2000 or Windows XP operating systems. It plugs into LabVIEW and MATLAB and provides an SDK for customizable application development.

3.3. Note on operating systems

The MotionPro Data Acquisition System is supported in the following operating systems: Windows XP (32-bit), Windows Vista (32/64-bit) and Windows 7 (32/64-bit).

3.4. Software Development Kit

Upon the installation of the MotionPro DAS SDK several options are available to the user. These options are easily accessed via the Program menu under the Windows Start button. The programs and associated files are organized under the **IDT/XsDA** folder. This folder includes the example programs and the associated documentation. The software components included in the MotionPro DAS Software Development Kit are:

- Data Acquisition stand-alone application.
- SDK modules with example source code in MSVC++.
- ActiveX Control.
- Plug-in for LabVIEW™.
- Plug-in for MATLAB™.

The SDK modules provide an API interface to develop applications to operate the device and access all the capabilities using a programming language such as C++ and Java. A C/C++ header file is included in the SDK (**XsdaAPI.h** file in the Include sub-directory).

For a more detailed description of the SDK please refer to the **MotionPro DAS SDK Reference**.

4. Installing the MotionPro DAS

This section specifies the minimum recommended computer requirements and gives the procedures needed to install the device, I/O Cable, and software.

4.1. Minimum computer requirements

	PC
Operating System	Windows XP, Vista, 7
Processor	Pentium III or equivalent with 500 MHz processor.
RAM	2 GB
USB 2.0 Port	high speed USB port that is NOT shared with other devices
Hard Drive	100GB or larger hard drive (recommended).

4.2. Package contents

Before beginning the installation process, check that the following items are present in the MotionPro DAS package. If you are missing any of the items listed below, please contact IDT, Inc. or your sales representative.

- Data Acquisition module.
- I/O USB 2.0 Cable.
- CD-ROM of Data Acquisition software and documentation.
- Data Acquisition Quick Start Guide.

4.3. Software Installation



Windows XP/Vista/7/8/8.1

Before installing the software make sure that the computer has Windows 2000 or Windows XP installed as operating system.

- Log into Windows with a username and password that has **ADMINISTRATIVE PRIVILEGES**.
- Insert the Data Acquisition CD in the CD drive. If the computer is configured to AUTORUN, the installer will run automatically. If not, click on the Windows **Start** button. Select **Run** from the menu. Use the **Browse** button to locate the SETUP.EXE file on the Data Acquisition CD and click the **OK** button.
- Select the "Install" option and follow the on-screen instructions.
- EXIT when the installation is complete and restart your computer.

4.4. Hardware Installation

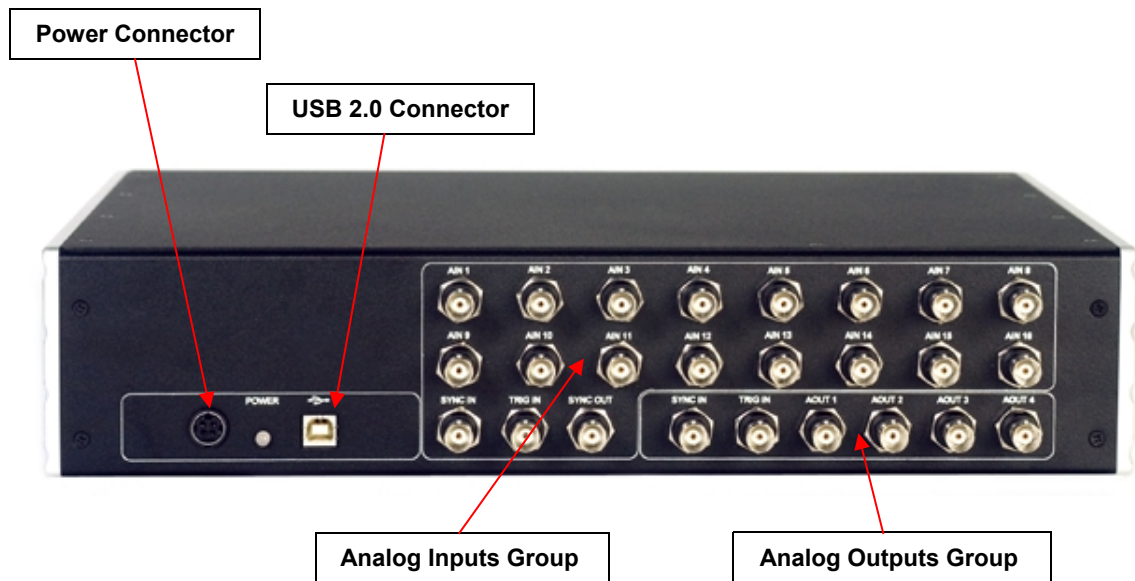
A 5 VDC, 6 Amp supplies the device with the necessary power. This power supply unit is included with the system package. All communication and data transfer with the host computer is done via the USB-2 interface. This interface requires a single cable, which is also supplied with the camera package.



Windows 2000/XP/Vista/7

- Connect the device to the power source.
- Connect the USB 2.0 cable to an available USB 2.0 port on your computer.
- Connect the other end of the USB 2.0 cable to the device and wait a few seconds for the device to initialize it self.
- Follow the on-screen instructions. Click on the **YES** or **Continue Anyway** button when prompted by the Windows 2000 or XP Operating system to proceed with the installation.

4.5. Data Acquisition Back panel



The BNC connectors are grouped in two groups:

Analog Inputs

- **AIN1 to AIN15:** the analog input channels
- **Sync In:** the analog input external synchronization input
- **Trig In:** the analog input external event trigger input
- **Sync Out:** the analog input output signal.

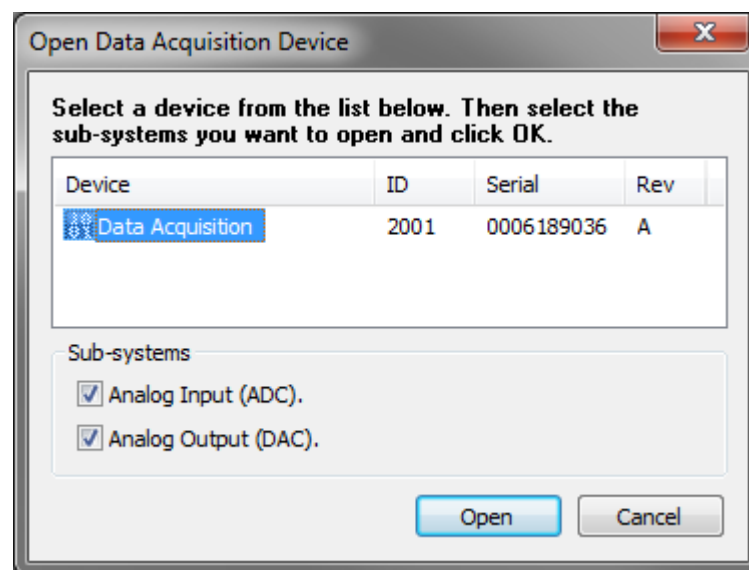
Analog Outputs

- **AOUT1 to AOUT4:** the analog output signals
- **Sync In:** the analog output external synchronization input
- **Trig In:** the analog output external event trigger input

5. Data Acquisition stand-alone program


This application allows the user to acquire data from the analog inputs, save the data to the hard disk, generate waveforms for the analog output channels, open saved data files and display the signals, apply filters and mathematical operations to the signals. Upon execution the user can select to open a data acquisition window.

5.1. Open data acquisition device



The “Open data acquisition device” dialog box allows the user to open and configure a data acquisition device. If more than one device is attached to the computer system, the user may select from a list. For each device, the user may select the analog subsystems that will be open.

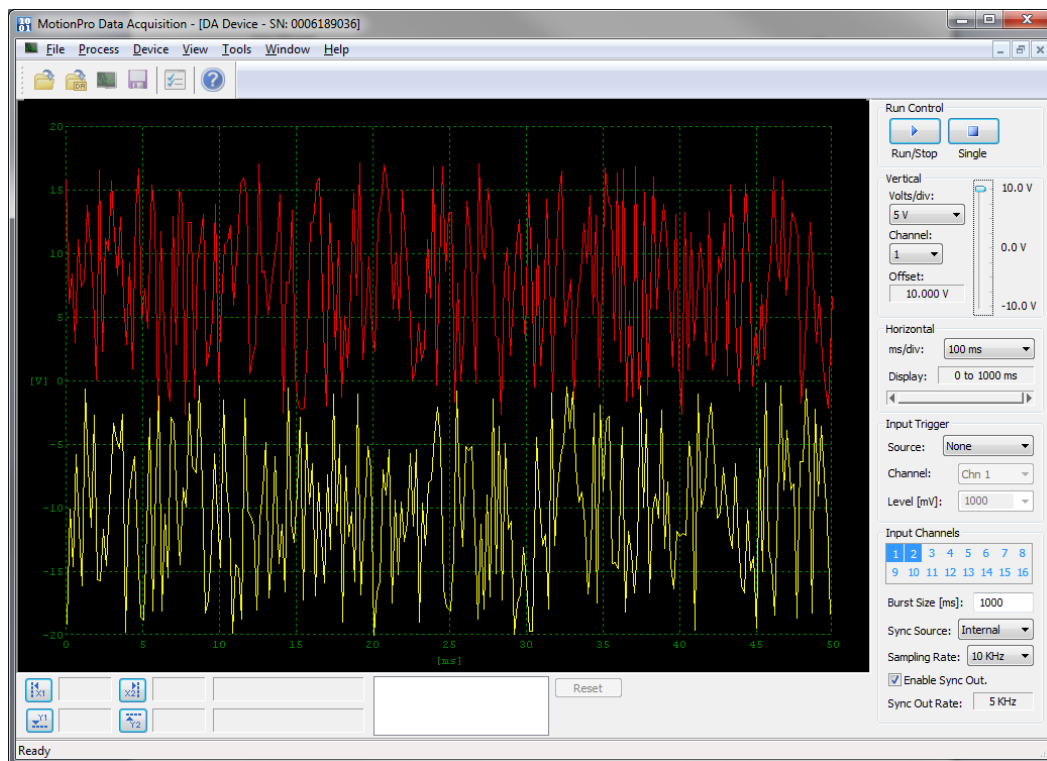Once the dialog box completes the initialization procedure for the device, the application window for the active device is open as shown below.



Along with the usual top bar menu structure the application also includes a docked dialog bar on the right side and a docked dialog bar on the bottom side. In the right bar the main operational controls of the data acquisition are grouped by function: Run control for acquisition/waveform generation; vertical and horizontal configuration for data display; input trigger for the configuration of the ADC trigger. The function of these controls is in great part also accessible from the top menu bar. In the bottom bar, four buttons allow the configuration of horizontal and vertical cursors for the computation of time, frequency and voltage. A list shows the filter and mathematical operations applied to the signals and a reset button allows the user to restore the original signal.

5.2. Data Acquisition Menu structure

The Data Acquisition main menu contains the following options:

File

Process

Device

View

Tools

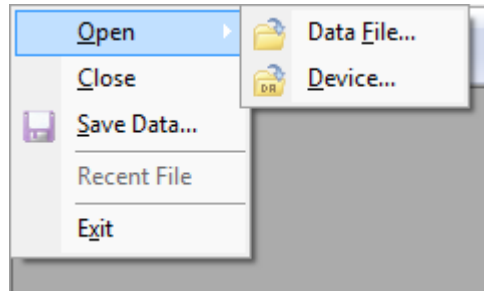
Window

Help

5.3. FILE Menu

The file menu contains the following options:


- Open the data acquisition window.
- Open a data file.
- Save acquired data on the hard disk and close windows.



5.3.1. Open a data file

Each stored data file may be open and displayed.

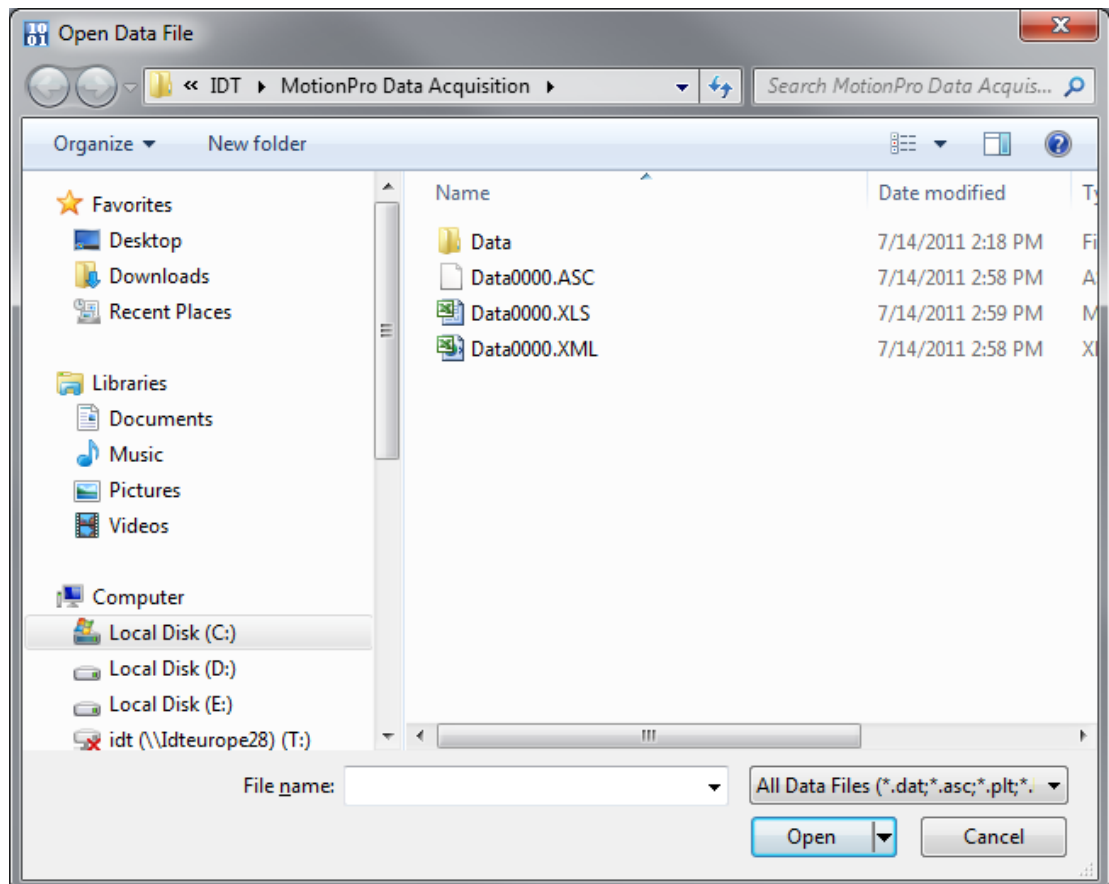1. From the File menu select File > Open > Data File...
2. Select the desired folder and the file name.
3. Click Open.



5.3.2. Save acquired data

Each acquired sequence of data may be stored on the hard disk.

1. From the File menu select File > Save Data.
2. Select the desired folder.
3. Type the file name, select the file format from the “Save as type” control.
4. Click “Save”.

The user may save data in the formats listed below:

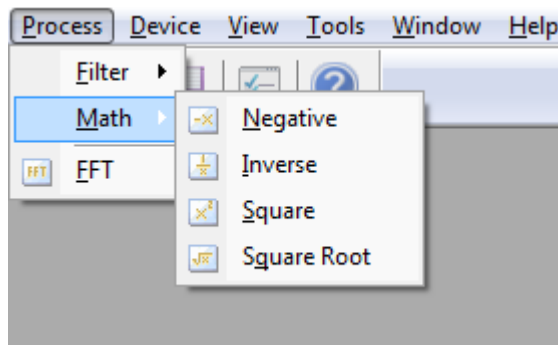
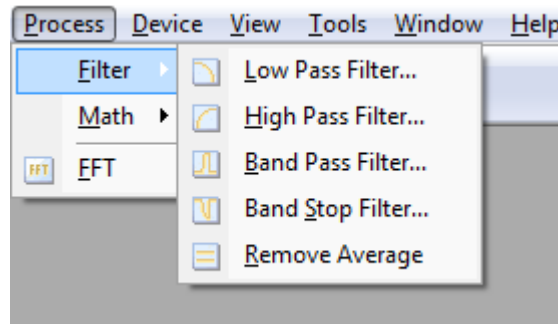
- ASCII File (TXT).
- Binary File (DAT).
- Tecplot File (PLT).
- LabVIEW Measurement File (LVM).
- Excel XML spreadsheet File (XML).

The values saved in a data file may be used as sources to generate output signals. For further information, refer to the output channels paragraph.

5.4. PROCESS Menu

The Process menu contains the following options:


- Apply filters and remove average from the signals.
- Apply mathematical operations to the signals.
- Show/Hide the FFT/Power Spectrum window.



5.4.1. Filters and Math operations

Filter operations: low pass filter, high pass filter, Band Pass filter, Band Stop filter or Remove average.

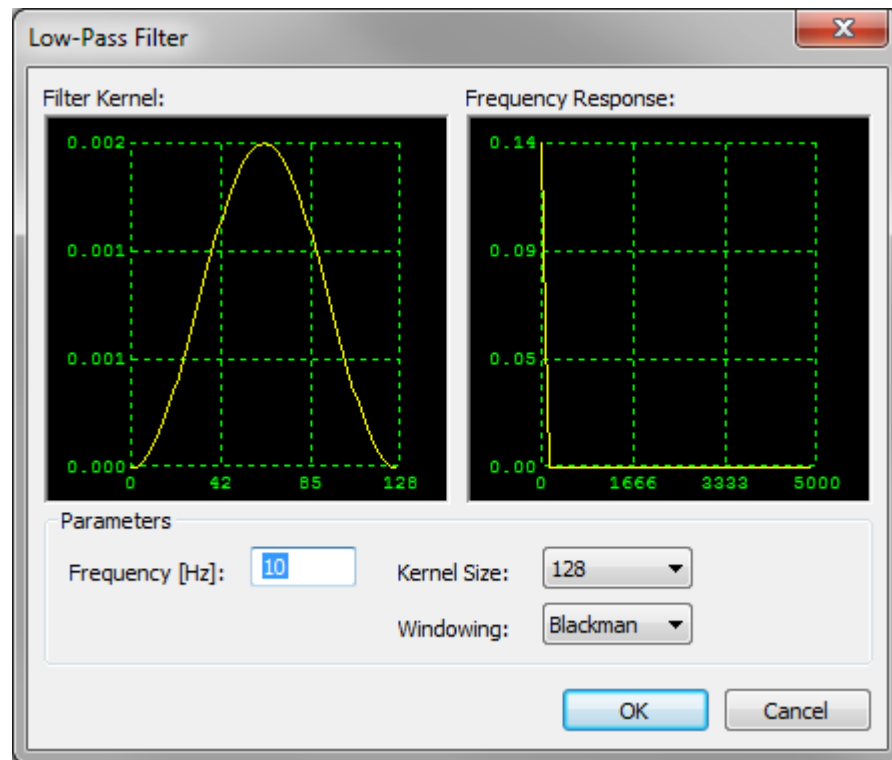Mathematical operations: Negative, Invert, Square and Square Root.



The filter dialog box shows the kernel and the Frequency response of the selected filter.

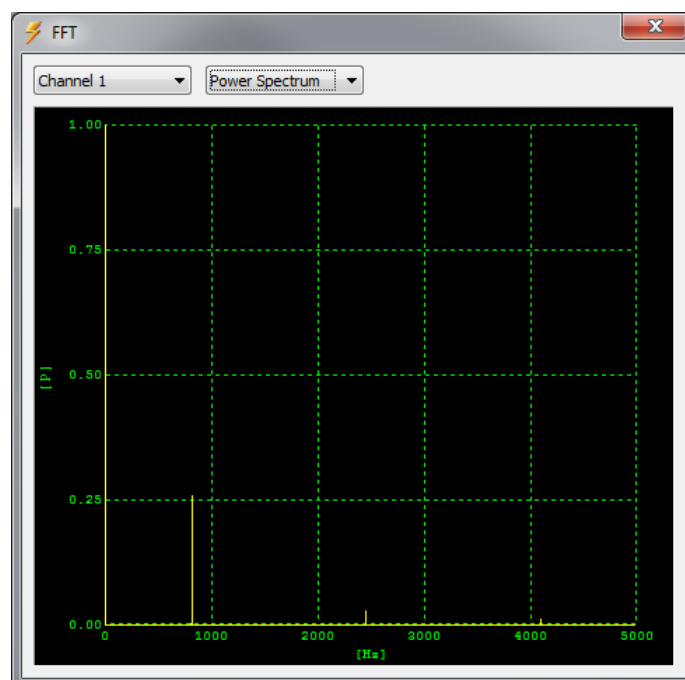
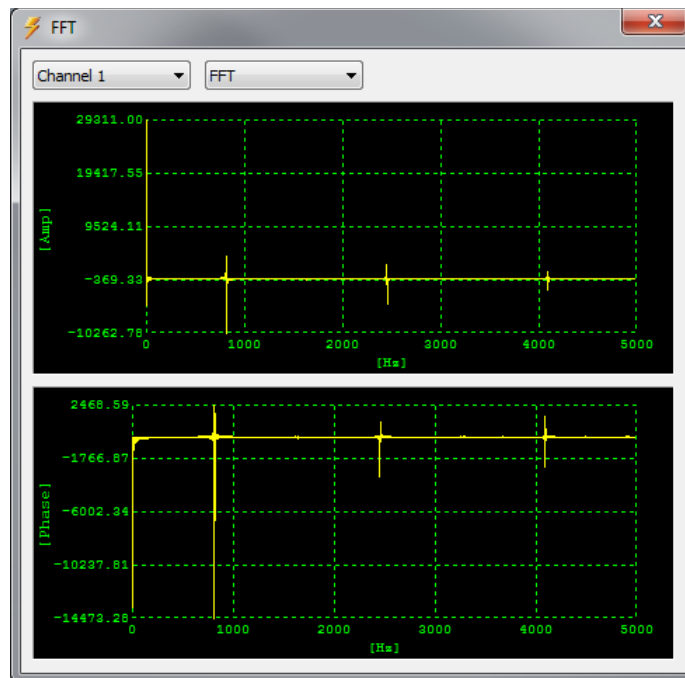
The user may select the frequencies (one value for low pass and high pass filters, two values for band pass and band stop filters), the kernel size and the windowing type (rectangle, Bartlett, Blackman, Hamming and Hann).

If the filters or the mathematical operations are activated the operations list in the horizontal bar is filled with those operations. The user may:

- Select the items in the list and activate the single operation.
- Reset the list and delete the operations.

5.4.2. FFT and Power Spectrum


If the “FFT” option is selected, the Fast Fourier Transform window is displayed. The user may select the channel to display, the FFT and the power spectrum.



5.5. DEVICE Menu

The Device Menu on the main toolbar offers an alternative to using the buttons provided by the Docked Dialog menu including the following functions:

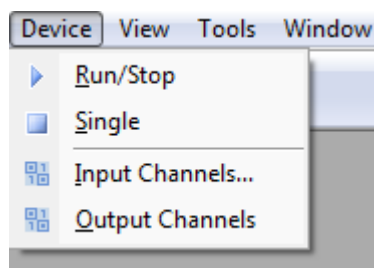1. Run/Stop.
2. Single.
3. Input channels.
4. Output channels.



5.5.1. Run Control

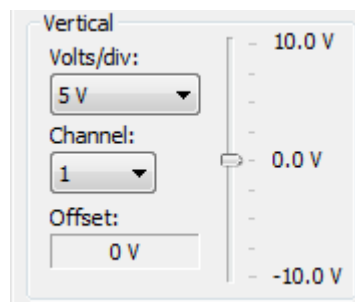
The Run Control group at the top of the Docked Dialog menu has data acquisition control functions including the following:

1. Run/Stop: start and stop continuous acquisition.
2. Single: run a single burst.



5.5.2. Vertical controls

Use the vertical controls group in the docked vertical bar to configure the data display vertical parameters.




Volts/div: the control sets the vertical resolution in volts/division.

Vertical displacement (Offset): the user may vertically scroll each input signal independently. If more than one signal is displayed, select the signal from the channel list and scroll the vertical slider.

5.5.3. Horizontal controls

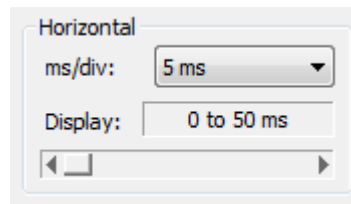Use the horizontal controls group in the docked vertical bar to configure the data display horizontal parameters.




ms/div: the control sets the horizontal resolution in ms/division. The option “Fit to window” automatically sets a value that fits the signal length in the current window.

Display: scroll the blue bar to display different parts of your acquisition. The bar is automatically resized if the length of the acquisition cannot fit in the screen.

5.5.4. Input trigger

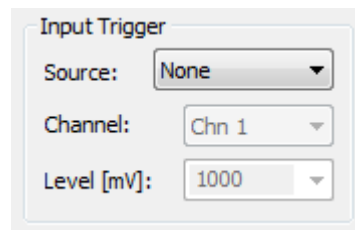The acquisition of data is started by a trigger that can be configured by the user.

A screenshot of a software dialog box titled "Input Trigger". It contains three rows of controls. The first row is "Source:" with a dropdown menu showing "None". The second row is "Channel:" with a dropdown menu showing "Chn 1". The third row is "Level [mV]:" with a dropdown menu showing "1000". The "Channel:" and "Level [mV]:" controls are grayed out, indicating they are inactive when the "Source" is set to "None".

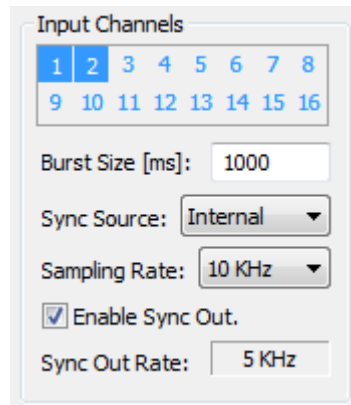
Source: the trigger source may be software (none), channel (one of the input channels themselves) or external (edge high or edge low).

Channel: if the trigger source is channel, use this control to select it. The option is grayed out if the source is not channel.

Level: use this option to select the trigger level value. The option is grayed out if the source is not channel.

5.5.5. Input channels

Use the input channels controls to configure the acquisition channels.



The screenshot shows the 'Input Channels' configuration window. It features a grid of 16 channel numbers (1-16) arranged in two rows. A blue rectangle highlights the first two channels (1 and 2). Below the grid, there are several controls: 'Burst Size [ms]' set to 1000, 'Sync Source' set to 'Internal', 'Sampling Rate' set to '10 KHz', a checked 'Enable Sync Out' checkbox, and 'Sync Out Rate' set to '5 KHz'.

Channels display: the rectangle shows the open channels highlighted by a blue rectangle. Click on the numbers to configure the channels.

Burst size: the size of the acquisition in ms.

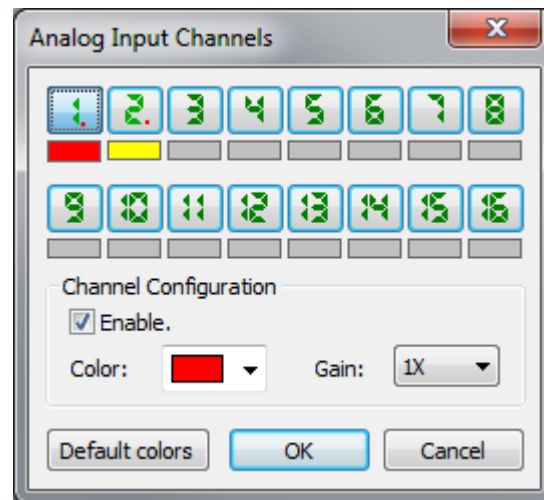
Sync source: the synchronization source (internal or external).

Sampling rate: select the sampling rate.

Enable Sync Out: check button to enable the sync out signal. The minimum number of input channels required to activate the sync out is 2. The sync out frequency is equal to the Sampling rate divided by the number of input channels.

5.5.6. Configuring input channels

If the input channels rectangle is clicked the dialog box below appears.



Channels button: click on the button of the channel that you want to configure.

Enable: select this option to enable/disable the channel. A disabled channel is not sampled.

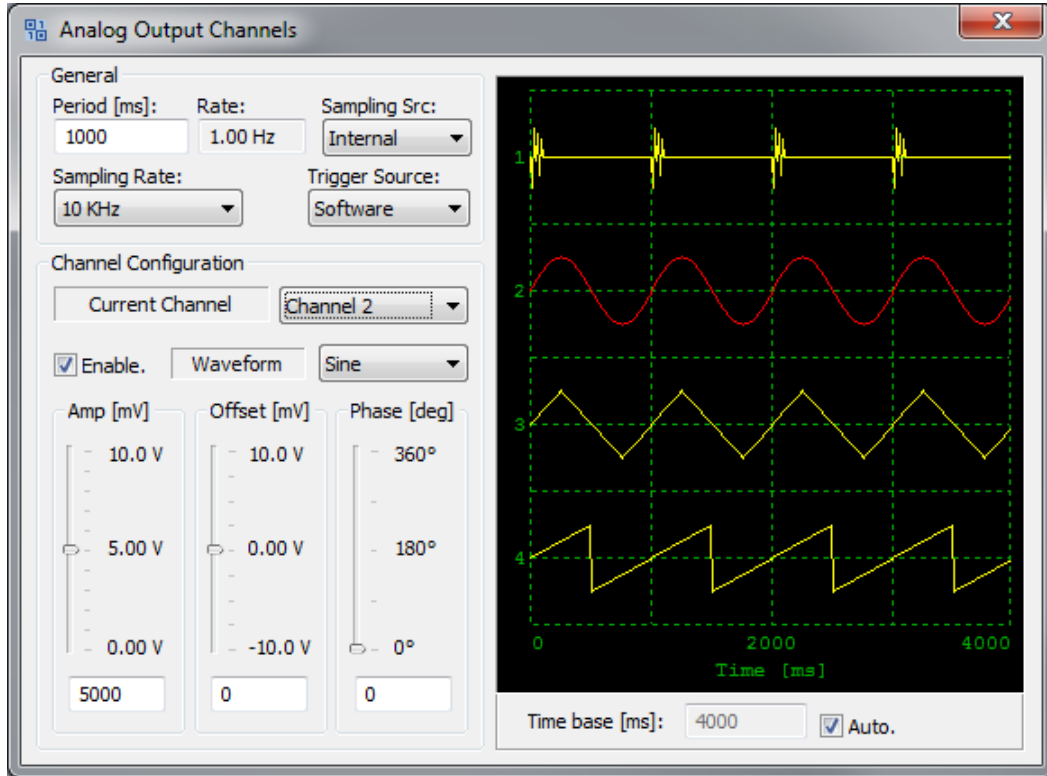
Color: select the channel display color.

Gain: select the channel gain (1X, 2X, 4X, 8X).

Default color: click this button to assign the default colors to the channels

5.5.7. Output channels

Click the output button on the docked vertical bar and the output channels dialog box will appear. Output parameters are independent from the input but they are activated by the play button on the vertical toolbar.



Period: the output signal period in ms.

Sync source: the output sync source (internal or external).

Sync rate: the output rate.

Trigger source: the start trigger source may be software or external (edge high or edge low).

Four analog output channels are available. Each channel may be independently configured.

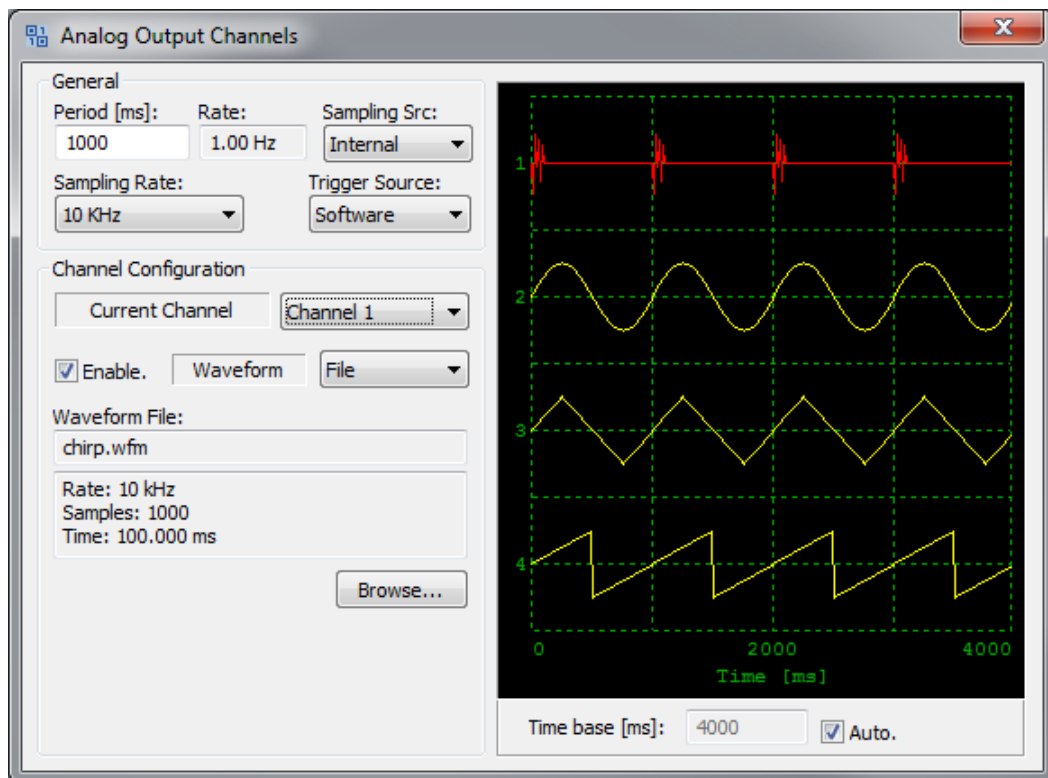
Current Channel: select the channel that will be configured.

Enable: select this option to enable/disable the channel. A disabled channel is not sent to the output DAC.

Waveform: select the waveform shape (square, sine, triangle, saw tooth or File)

Amplitude - Offset - Phase: configure the waveform parameters (square, sine, triangle, saw tooth).


File: if the waveform type is file the output channel generates a waveform read from a WFM file. WFM files are generated and edited by the Waveform Editor. For further information, please refer to the chapter 6.



The signals are shown in the black window. The user may select the Time Base period in μ s. If the "Auto" option is selected, the application shows four waveform periods.

5.6. View Menu

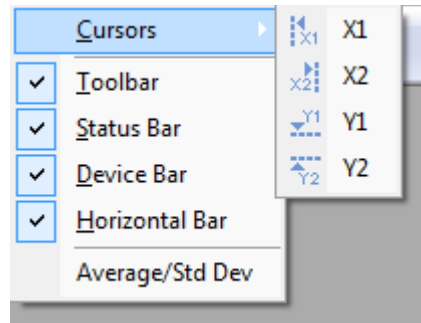Use the View menu to activate the vertical (X1 and X2) and the horizontal (Y1 and Y2) cursors. Use the View menu to select the toolbar options.



5.6.1. Average and Standard Deviation

The "Average/Std Dev" option in the View menu activates the "Average/Standard Deviation" window.

A screenshot of the 'Average / Standard Deviation' window. The window has a title bar with the text 'Average / Standard Deviation' and a red close button. Inside the window is a table with five columns: 'Channel', 'Minimum', 'Maximum', 'Average', and 'Standard Dev'. There is one row of data for 'Channel 1'.


Channel	Minimum	Maximum	Average	Standard Dev
1	500.3 mV	5479.7 mV	2951.0 mV	2413.4 mV

For each channel the minimum, the maximum, the average and the standard deviation values are displayed.


5.6.2. Cursors

You may measure waveform data using cursors. Cursors are horizontal and vertical markers that indicate X-axis values (time) and Y-axis values (voltage) on a selected waveform source.

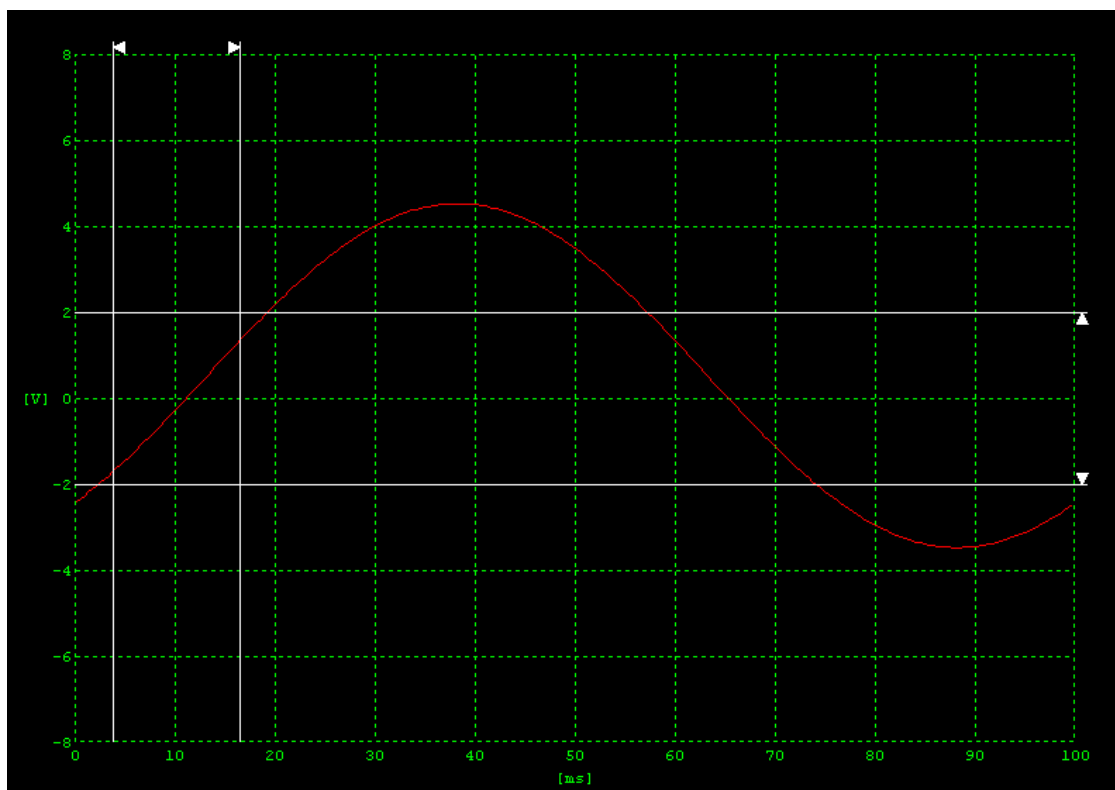The position of the cursors can be moved by dragging the triangular handles on the screen. See the picture below.



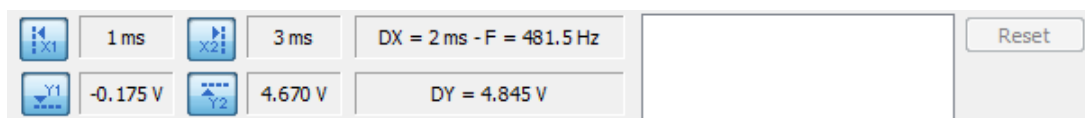The cursors are activated by clicking the corresponding buttons in the horizontal bottom bar.



In the same bar the displayed values are $DX (X_2 - X_1)$, the frequency and $DY (Y_2 - Y_1)$.

Cursors are not always limited to the visible display. If you set a cursor, then pan and zoom the waveform until the cursor is off-screen, its value will not be changed, and if you pan the waveform back again it will have the cursor in its original position.

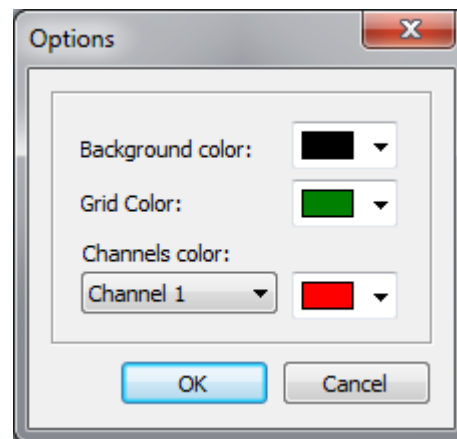
5.7. TOOLS menu

5.7.1. Waveform Editor

1. From the Tools menu, select “Waveform Editor...”
2. The waveform editor is launched; please refer to chapter 6.


5.7.2. Options

1. From the Tools menu, select “Options...”
2. Select the colors and click OK.



5.7.3. Language selection

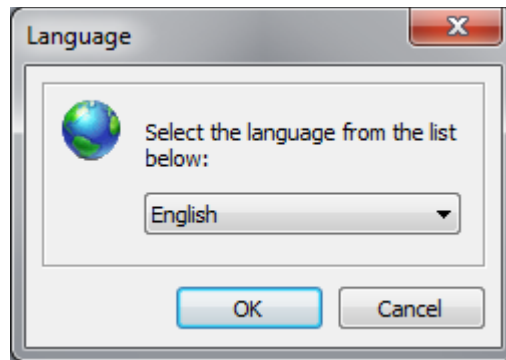1. From the Tools menu, select “Language...”
2. Select the language from the list and click OK.



5.8. WINDOW Menu

If more than a window is open in the program, use the window menu to cascade, tile horizontally, tile vertically or select one of the windows.

5.9. HELP Menu

This menu contains support options and information including: e-mail tech support and software and manual updates.

6. Waveform Editor

6.1. Overview

The waveform editor is a tool for the generation of analog waveforms files.

A **waveform** is made up of one or more components.

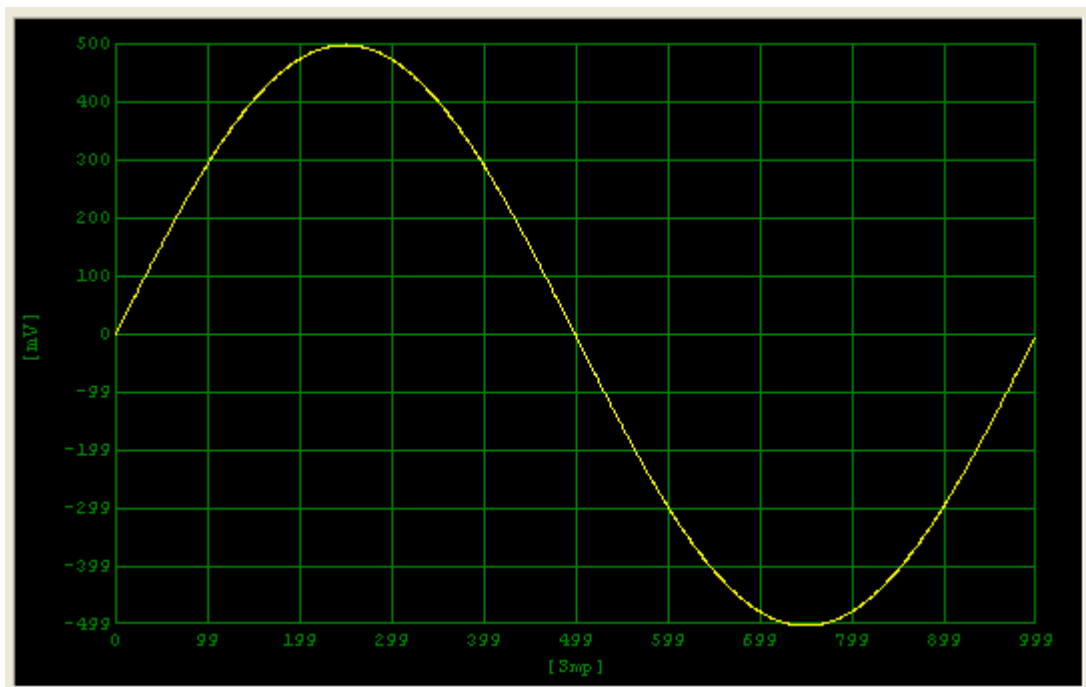
A **component** is made up of one or more primitives.

A **primitive** is the lowest level of a waveform. A primitive can be anything from sine waves and triangle waves to pulses and noise. Primitives of the component are added, subtracted, multiplied, or divided into each other.

The examples below show different waveforms.


Example #1

The example below shows a waveform with one component. The component has a single sine wave primitive.



Example #2

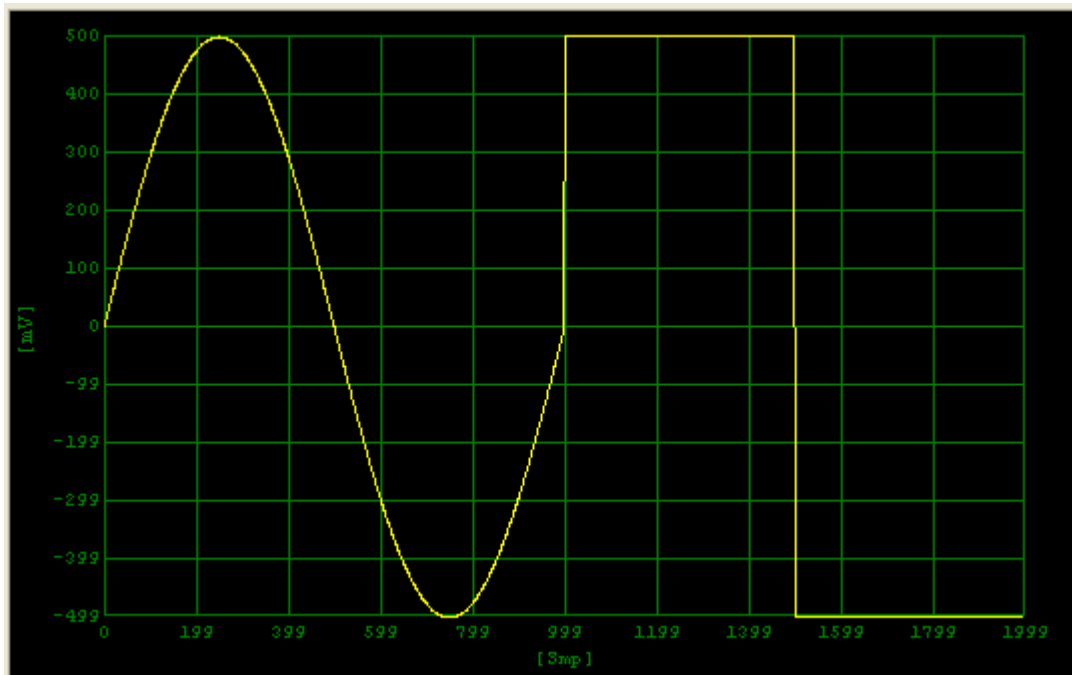The example below shows a waveform with two components. Each component is made up of one primitive, a sine and a square wave.



Example #3

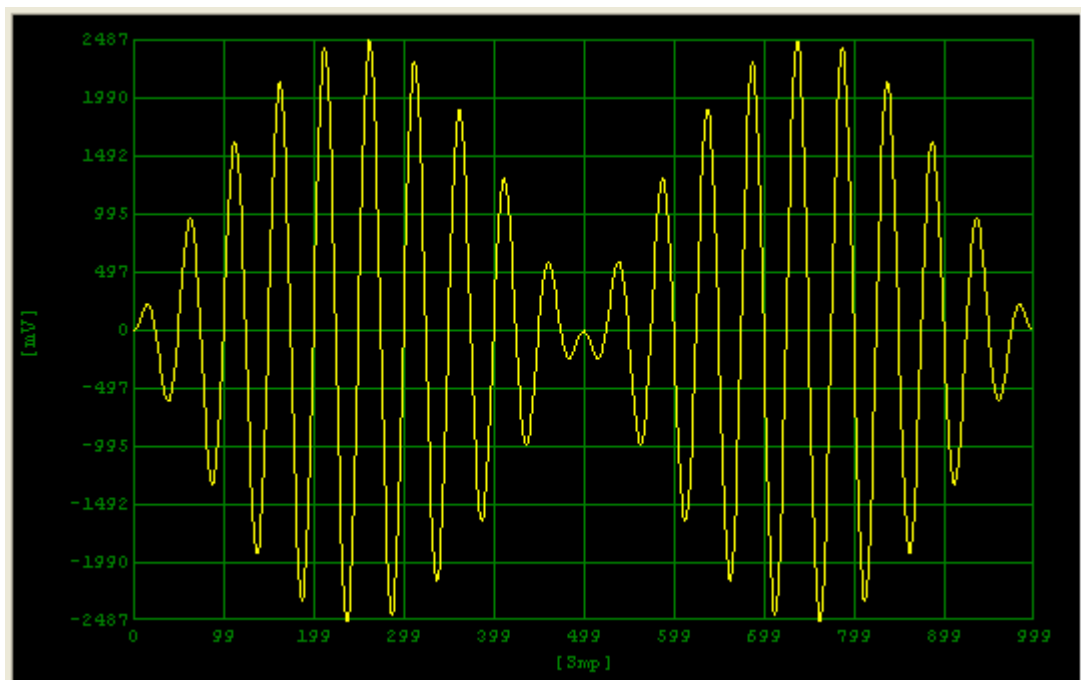
The example below shows a waveform made up of one component. The component is made up of 2 primitives.

The **first primitive** is a sine wave with the following properties:

- Amplitude: 500 mV
- Period: 1000 Samples
- Phase: 0°

The second primitive is a sine wave with an associated operation of multiplication. This primitive has the following properties:

- Amplitude: 5 mV
- Period: 50 Samples
- Phase: 0°



6.2. The waveform editor menu structure

The Data Acquisition main menu contains the following options:

File

Edit

View

Tools

Help

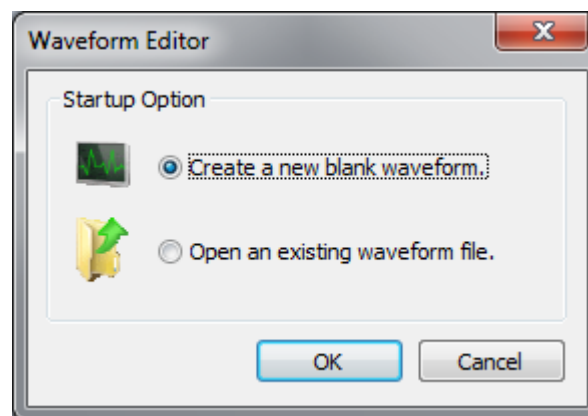
6.3. File Menu

The file menu contains the following options:


- Create a new waveform.
- Open a waveform file.
- Save a waveform file.

6.4. Create a new waveform

Upon execution of the program the user can select to create a new waveform or open an existing waveform file.



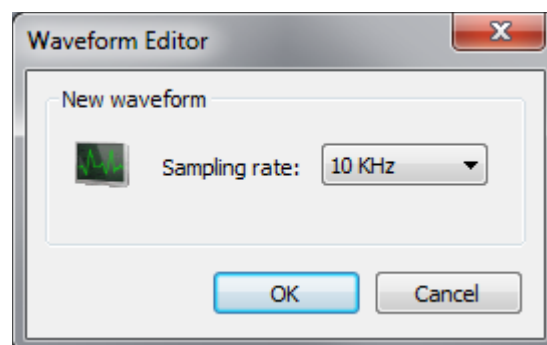1. Select the "Create a new blank waveform" option and click OK.
2. Select the sampling rate and click OK.




6.5. Open a waveform file

1. Upon execution of the program, select the “Open an existing waveform file” option and click OK.
2. Browse for the waveform file to open.
3. Click the Open button.

6.6. Save a waveform

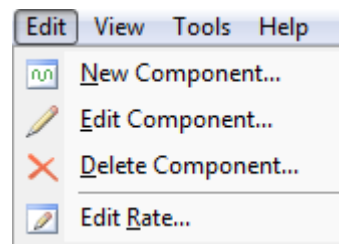
 To save a waveform file, from the File menu, select the “Save waveform...” item.

 Type the file name and click OK.

6.7. Edit Menu


The edit menu contains the following options:

- Add a new component to the waveform.
- Edit one of the components of the waveform.
- Delete one of the components of the waveform.
- Edit the waveform sampling rate.



6.8. Add components to a waveform

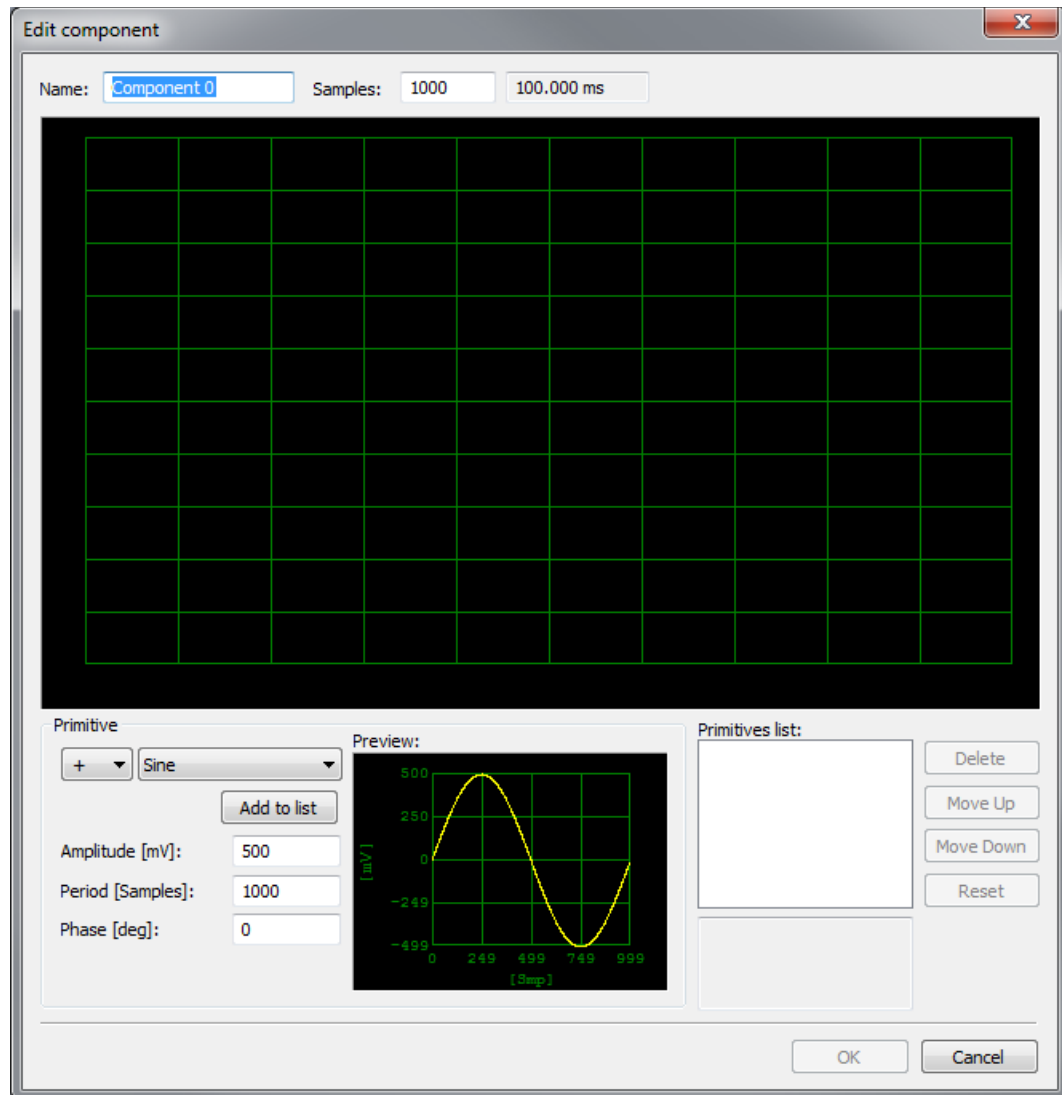
1. From the “Edit menu” select the “New Component...” item.
2. The component editor window opens.



You may select the component name and the number of samples. The samples are converted into ms according to the selected sampling rate.


6.9. Insert primitives into a component

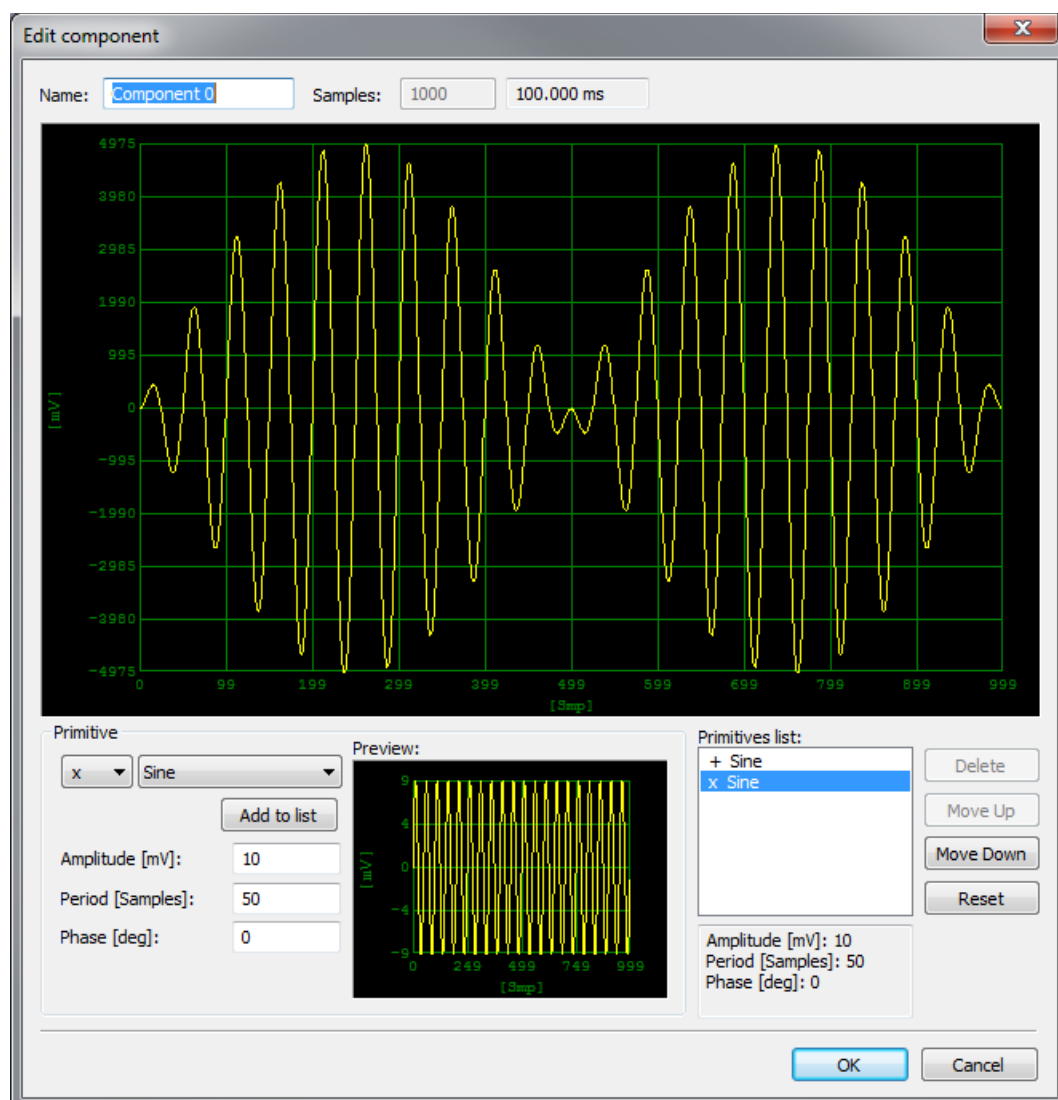1. From the component editor window select the primitive type and operation (+, -, x, /).
2. Modify the primitive parameters. Different primitive types have different parameters.
3. Click the “Add to List” button.
4. You may also:
 - Delete a primitive from the list.
 - Reset the list
 - Move a primitive Up and Down in the list.




6.10. Edit a component

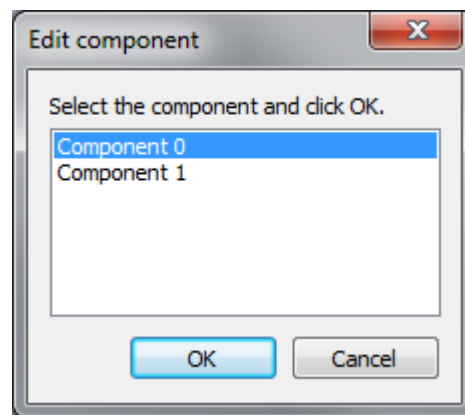1. From the Edit menu select the “Edit Component...” item.
2. If the waveform has a single component the edit component window will appear, otherwise a new window will ask you to select which component to edit.
3. Select the component from the list and click OK.
4. Refer to the previous topic to add, delete and edit the primitives in the component.



6.11. Delete a component

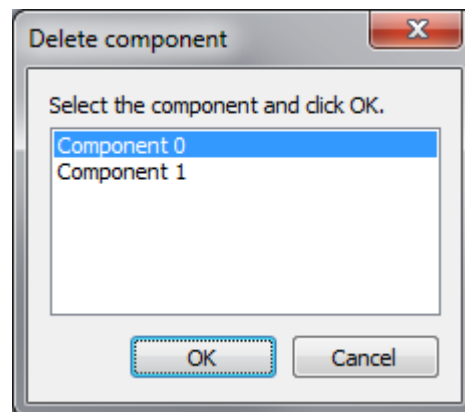1. From the Edit menu select the “Delete component” item.
2. If the waveform has a single component the component will be deleted, otherwise a new window will ask you to select which component to delete.
3. Select the component from the list and click OK.

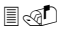



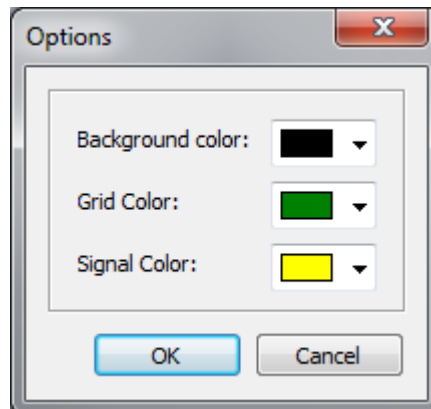
6.12. View Menu

Use the View menu to show or hide the main toolbar and the status bar.

6.13. Tools Menu

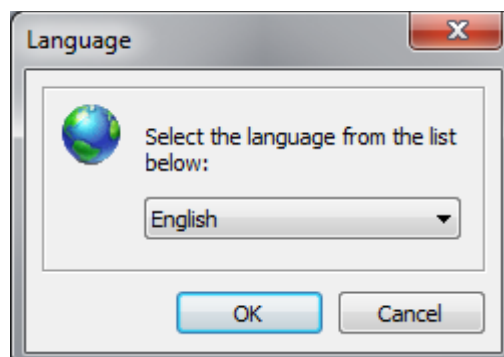
6.13.1. Options

-  From the “Tools” menu, select the “Options...” item.
-  Change the background, grid, and signal colors, then click OK.



6.13.2. Language

1. From the Tools menu, select “Language...”
2. Select the language from the list and click OK.



6.14. Help Menu

This menu contains support options and information including: e-mail tech support and software and manual updates.

7. Appendix A – Product Specifications

7.1. Specifications

Analog inputs	
Number of inputs	16
Resolution	16 bit
Programmable gains	1X, 2X, 4X, 8X
Ranges	$\pm 10\text{ V}$, $\pm 5\text{ V}$, $\pm 2.5\text{ V}$, $\pm 1.25\text{ V}$
Single Channel A/D throughput	500 kS/s
Multi channel A/D throughput	500 kS/s $\pm 0.05\%$
A/D conversion time	2 μs
Bipolar input range	$\pm 10\text{ V}$
Output coding	Offset binary
Maximum input voltage (power off)	$\pm 35\text{ V}$
Maximum input voltage (power on)	$\pm 20\text{ V}$
Input impedance (channel off)	100 M Ω , 10 pF
Input impedance (channel on)	100 M Ω , 100 pF
Bias current	$\pm 20\text{ nA}$
Non linearity	$< 1/2\text{ LSB}$
Analog outputs	
Number of DACs	4
Resolution	16 bit
Throughput	500 kS/s
Output range	$\pm 10\text{ V}$
Data coding	Offset binary
Output current	$\pm 5\text{ mA}$
Output impedance	0.1 Ω
Capacitive driver capability	0.004 μF
Protection against	Short circuit to analog ground
Non linearity	1 LSB
External A/D and D/A Triggers	
Internal trigger source	Software initiated
External trigger source	Software selected
Input type	Edge sensitive
Logic	TTL
Input termination	2.2 k Ω pull-up to +3.3 V
Logic high input voltage	2.0 V
Logic low input voltage	0.8 V

Minimum pulse width	25 ns
Miscellaneous	
Interface	USB 2.0
Operating systems	Windows 2000/XP
Plug' n 'play	Yes
Power requirements	5 V DC at 2 A maximum
Operating temperature range	0 ° to +55° C
Relative humidity	95%